

SimpliShade™

FAMILY



SimpliShade™

UNIVERSAL

Simplified Universal Composite System with just three shades

SimpliShade covers all 16 VITA® classical shades with just 3 shades. Simplify your inventory management and save yourself time on shade selection. Keep your patients confidently smiling while keeping your workflow simple.



SimpliShade™

BULK FILL

Simplify with One Shade and One Step

Beautiful restorations just got simpler, with SimpliShade™ Bulk Fill. You can skip the capping layer with a one-step placement and reduce inventory requirements with one shade for all 16 VITA® classical shades.



SimpliShade™

BULK FILL FLOW

Single-Shade Simplicity Meets Superior Adaptability

Reach every area of the tooth's anatomy while saving time with SimpliShade Bulk Fill Flow. Requiring just one shade to match all 16 VITA® classical shades, SimpliShade Bulk Fill Flow easily adapts to the cavity floor — with exceptional durability and depth of cure.

SimpliShade™ UNIVERSAL

WITH JUST THREE SHADES

Use SimpliShade with confidence in all cavity classes without the need for a blocker.



Simplify Shade-matching

Shade selection made easier and faster. SimpliShade matches all 16 VITA classical shades with just 3 shades: Light, Medium, Dark.



Streamline Inventory Management

Save time and money with minimal shades to meet almost all your composite needs. Declutter your shelves, avoid the burden of re-ordering too many different part numbers, and reduce or eliminate expired products.






Adaptive Response Technology (ART)

The ART filler network, a proven technology found in Harmonize®, helps you achieve lifelike restorations with exceptional strength and unmatched esthetics.



Universal For All Cavity Classes

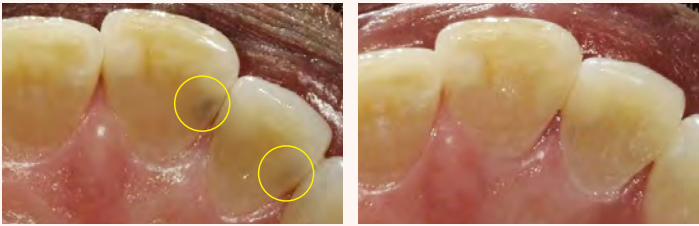Use SimpliShade with confidence in all cavity classes without the need for a blocker.


SimpliShade Before & After		VITA Shade Guide Equivalent
LIGHT MEDIUM DARK		 A1, B1, B2, C1, D2
		 A2, A3, C2, D3, D4
		 A3,5, A4, B3, B4, C3, C4

Case images courtesy of Dr. Matthew Miller



ART³: A NEW WAY TO BULK FILL

Achieve lifelike restorations with the next generation of clinically proven Adaptive Response Technology (ART) found in Harmonize™ and OptiShade Universal Composite. Powered by Enhanced Adaptive Response Technology, SimpliShade Bulk Fill offers easy handling and low shrinkage.

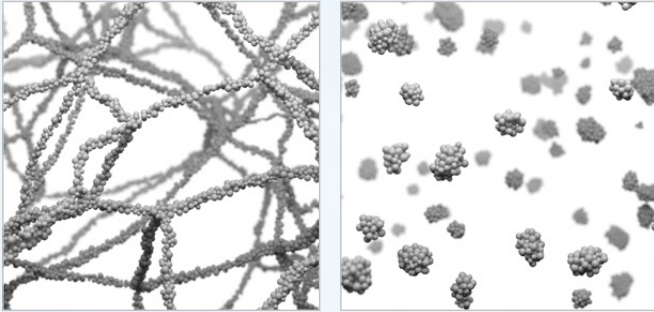


Adaptive Rheological Response Technology for Ease of Sculpting

ART³ uses rheological modifiers that reduce the composite's viscosity as you manipulate it. When you're done shaping, viscosity recovers and anatomy stabilizes before light-curing — reducing the need for resculpting.


Adaptive Color-Blending Response Technology for Beautiful Blending

Achieve exceptional blending with special nanofillers that mimic enamel structures and interact with light in a natural way. The restoration and tooth become virtually indistinguishable.

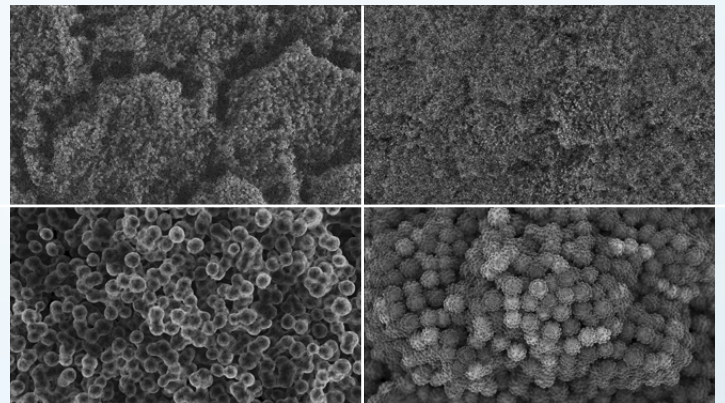

Human Enamel

SimpliShade Bulk Fill

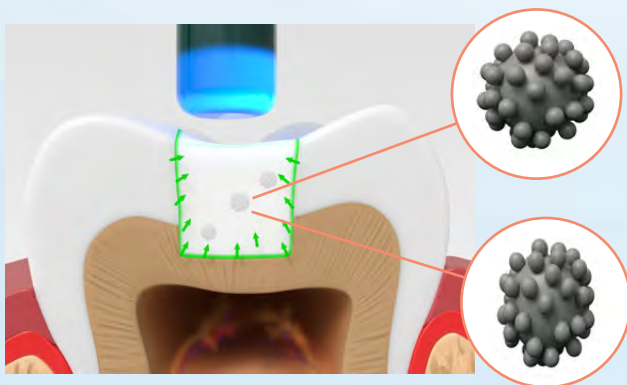

Adaptive Shrinkage Stress Response Technology for Durable Results

During polymerization, SimpliShade Bulk Fill uses low shrinkage resin and innovative nanofiller design, resin-infused nanofillers boost elasticity for stress relief and a special coating eases interfacial stress. Uniting three shrinkage reduction mechanisms, this enhances marginal adaptation and integrity by minimizing composite restoration stress on the adhesive layer for enhanced marginal adaptation.



Learn more about ART³ and the effect on shrinkage stress

SimpliShade™ BULK FILL

ONE SHADE IS ALL IT TAKES

SimpliShade™ Bulk Fill is indicated for Class I-V restorations - direct placement of anterior and posterior restorations, core build-ups, and repair of defects in restorations.



One-Shade System

Using one shade that blends seamlessly with all 16 VITA® classical shades, SimpliShade Bulk Fill saves you time selecting shades and managing inventory.



One Step

Eliminate the need for a capping layer, complete restorations with a SingleFill™ composite. With its enhanced curing mechanism, SimpliShade Bulk Fill makes it easy to achieve polymerization depths up to 5mm in a single step.



Simplicity with Durability

SimpliShade Bulk Fill offers low wear rate, high flexural strength, and comparable shrinkage stress relative to other bulk fill brands on the market, keeping your restorations stronger and long-lasting.





Chameleon Effect

Special nanofiller architecture mimics the natural enamel structure and interacts with light for virtually indistinguishable blending.

Flexural Strength


SimpliShade Bulk Fill shows slightly lower flexural strength compared to Filtek™ One Bulk Fill but higher flexural strength than Tetric® PowerFill and Admira® Fusion x-tra, resulting in better long-term resistance to chipping and higher fracture resistance.

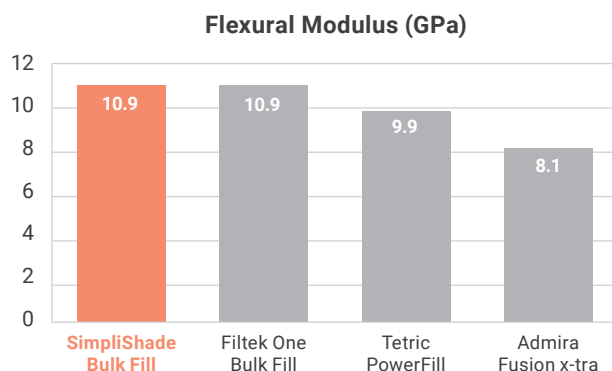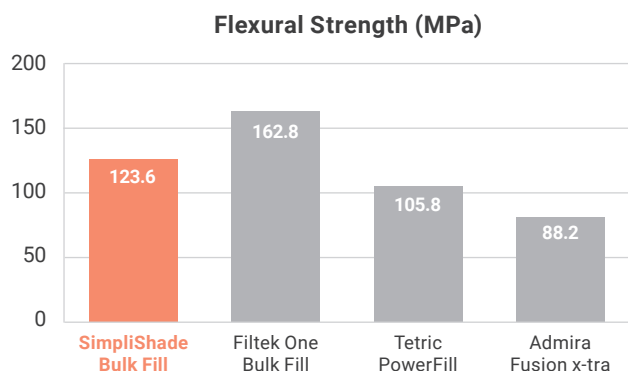

Authors: Drs. Jack L. Ferracane and Carmem S. Pfeifer, Oregon Health & Science University, Portland, OR.

Gloss Retention

SimpliShade Bulk Fill has the highest gloss retention compared to Filtek One Bulk Fill and Admira Fusion x-tra. This means that the esthetic initial gloss of the restoration is maintained for longer.

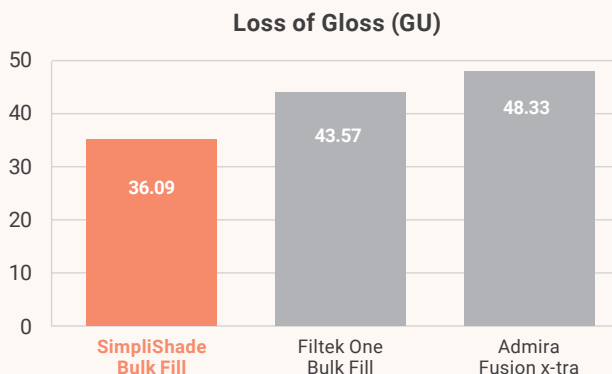

Authors: Prof. Dr. Nate Lawson, University of Alabama at Birmingham, Birmingham, AL.

Customer Testimonials



"SimpliShade Bulk Fill revolutionizes my approach to composites in molar and premolar restorations. Its convenience and ease of use with basic dental instruments are unmatched."

Dr. Maciej Mikołajczyk*
Lodz, Poland



"SimpliShade Bulk Fill offers a three-in-one solution with its universal shade, ideal consistency and excellent adaptability for simplicity and efficiency."

Prof. Joseph Sabbagh*
Brussels, Belgium



"Composite with wonderful consistency that allows easy application and modeling. Quick and simple work protocol, great shading."

Dr. Dimitar Kosturkov*
Sofia, Bulgaria

* Paid consultant for Kerr

SimpliShade™ BULK FILL FLOW

ONE SHADE IS ALL IT TAKES

SimpliShade™ Bulk Fill Flow is indicated for base/liner under **Class I** and **Class II** direct restorations, pit and fissure sealant, restoration of small cavities, **Class III** and **Class V** restorations, blocking out undercuts, small core build-ups, and repair of resin and acrylic provisional restorations.



Flowable Features Single-Shade Simplicity



One-Shade System

Using one shade that blends seamlessly to the cavity walls with all 16 VITA® classical shades, SimpliShade Bulk Fill Flow saves you time selecting shades and managing inventory.



Save Time

Simplify your workflow with just one step to cure up to 5 mm, before adding a capping layer.



Go with the Flow

Fill every area of the tooth's anatomy with a self-leveling bulk fill flowable that easily adapts to the cavity floor.



Low Shrinkage

The low polymerization-related shrinkage and shrinkage stress of SimpliShade Bulk fill Flow ensure long-lasting restorative results.




Add Confidence

Restore confidently with a flowable bulk fill composite that offers comparable or higher mechanical strength than the competitors.

Shrinkage Stress and Polymerization Stress


In terms of shrinkage behavior, there were no significant differences between the materials tested. The two competitor flowable materials, SureFil SDR Flow+ and Filtek Bulk Fill Flowable Restorative, are flowable composites with low shrinkage, so a similar performance of the SimpliShade Bulk Fill Flow material is a desirable result.

Shrinkage Stress (N)

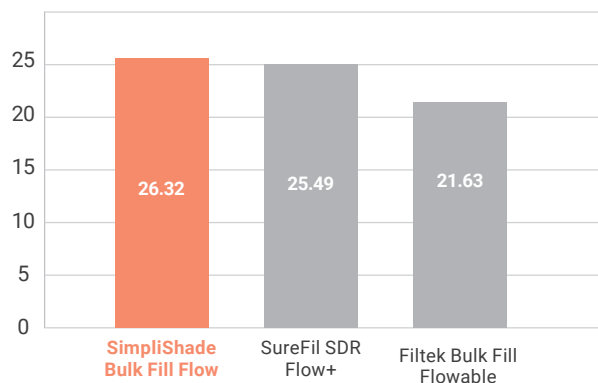


Authors: Prof. Dr. Nate Lawson, University of Alabama at Birmingham, Birmingham, AL

Polymerization Stress (MPa)

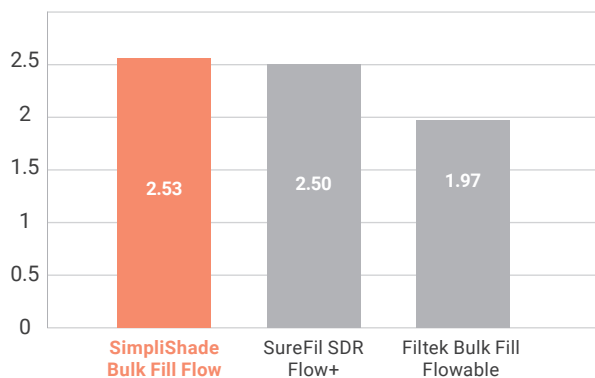

Authors: Drs. Jack L. Ferracane and Carmem S. Pfeifer, Oregon Health & Science University, Portland, OR

Volumetric wear (mm³)

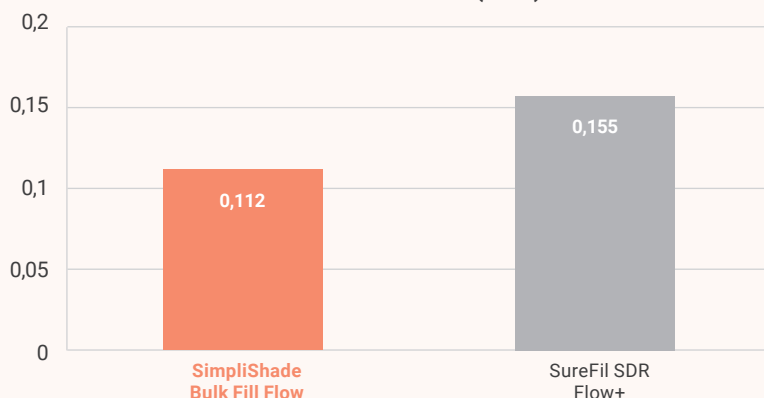

Authors: Prof. Dr. Nate Lawson, University of Alabama at Birmingham, Birmingham, AL

Wear Rate

SimpliShade Bulk Fill Flow offers a lower wear rate on the interproximal contacts than SureFil SDR Flow+, resulting in a long-lasting restoration.

Customer Testimonials



"Exceptional handling perfectly complements the Snow-plow restoration technique from Styleitaliano."

Prof. Angelo Putignano*
Ancona, Italy



"An impressive and innovative bulk fill flowable composite with great marginal adaptation, fast application with self-leveling and low shrinkage characteristics."

Dr. Iulian-Costin Lupu*
Iași, Romania




"Exceptional combination of good handling, low translucency and voids free flowable bulk fill composite with great radio-opacity."

Dr. Hossam Kamel*
Alexandria, Egypt

* Paid consultant for Kerr



Using SimpliShade Bulk Fill Flow (1), top up with SimpliShade Bulk Fill up to 5 mm (2) or with SimpliShade Universal.

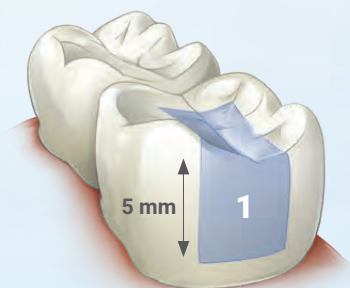

Using only SimpliShade Bulk Fill – Simply and effortlessly fill cavities up to 5 mm. No capping layer needed.

KEEP IT IN THE FAMILY

with our full line of SimpliShade products!



SimpliShade Universal Composite

Unidose

Part #	Item Description
37006	Light - 20 Pack
37007	Medium - 20 Pack
37008	Dark - 10 Pack
37009	Bleach White - 10 Pack
37010	Universal Opaque - 10 Pack

Syringe

Part #	Item Description
37001	Light
37002	Medium
37003	Dark
37004	Bleach White
37005	Universal Opaque



SimpliShade Bulk Fill

Part # Item Description

37211	SimpliShade Bulk Fill Composite Syringe - 4g
37212	SimpliShade Bulk Fill Composite Unidose 20 Pack x 0.25g




SimpliShade Bulk Fill Flow

Part # Item Description

37217	SimpliShade Bulk Fill Flow Syringe - 2x 2g
37218	SimpliShade Bulk Fill Flow Unidose 20 Pack x 0.25g



For an ultra high shine finish, use Kerr's line of Polishers!

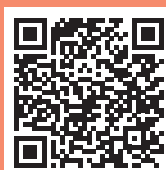

**One Shade.
One Step.
A Better Way
to Bulk Fill.**

For more information scan the QR Code or visit to.kerrdental.com/en/simplishadebulkfill